


# Köhler Illumination



# Köhler illumination

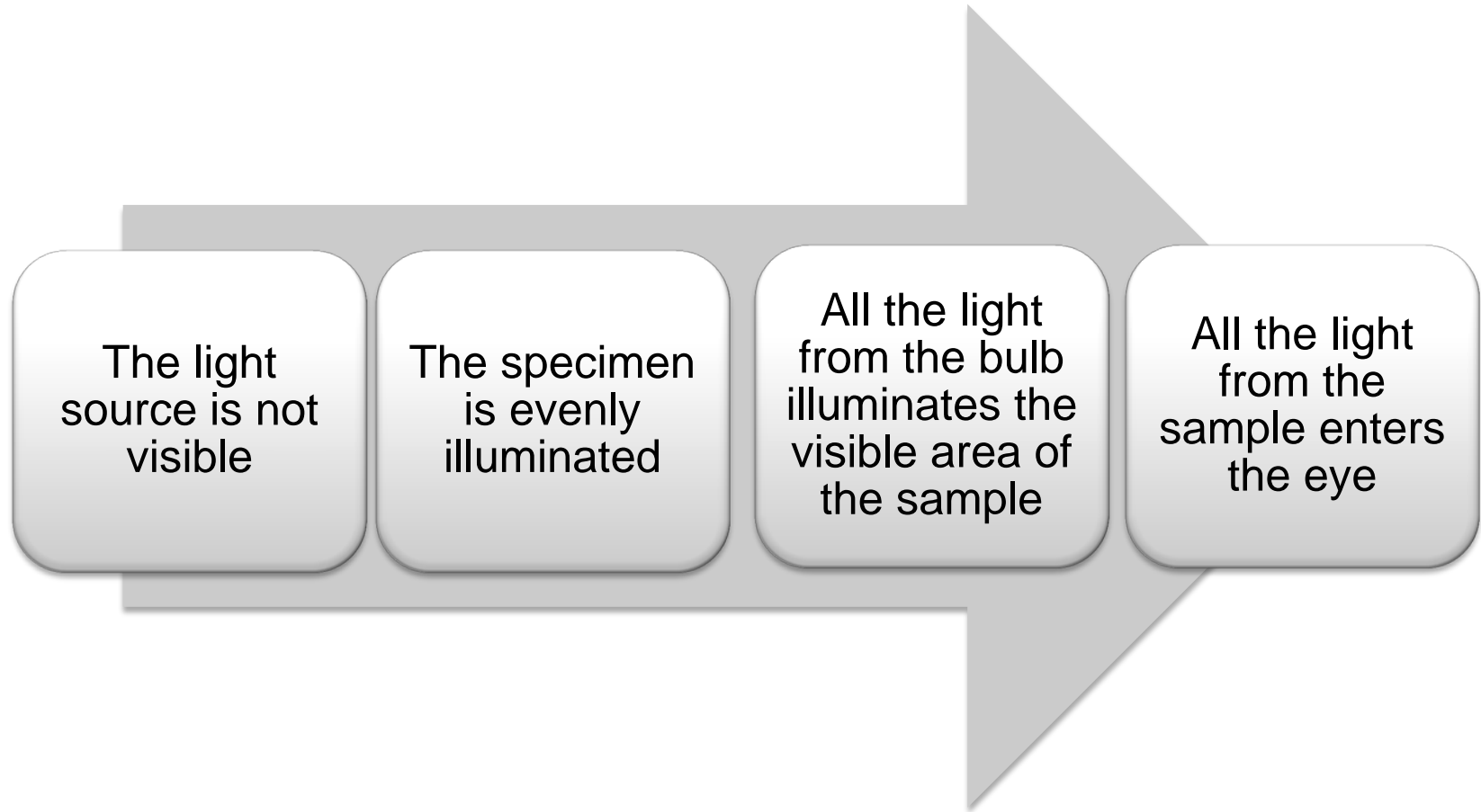

# Köhler illumination



Select 10X objective



Focus the sample

Close field diaphragm

Focus condenser

Open and centre the field diaphragm

Optimise condenser aperture

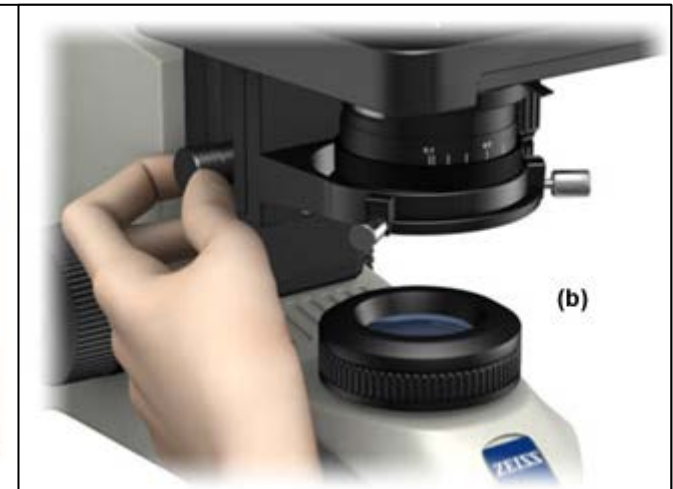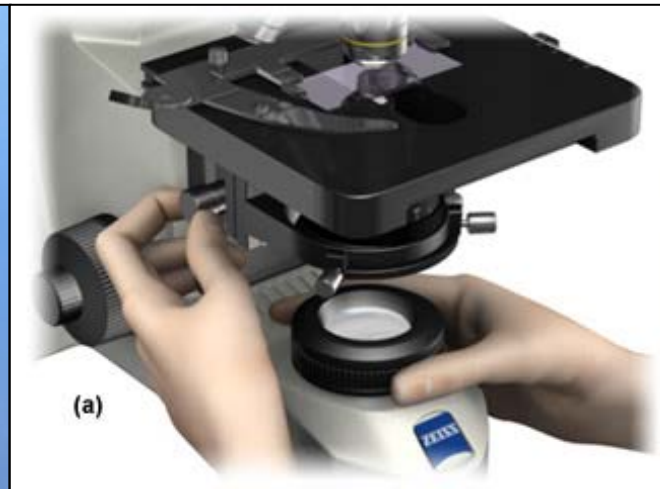

# Köhler illumination



Select 10X objective


Focus the sample

Close field diaphragm

Focus condenser

Open and centre the field diaphragm

Optimise condenser aperture

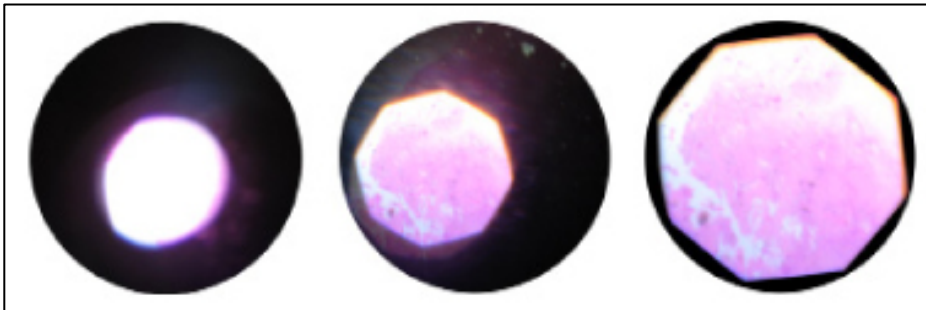

# Köhler illumination



Select 10X objective


Focus the sample

Close field diaphragm

Focus condenser

Open and centre the field diaphragm

Optimise condenser aperture

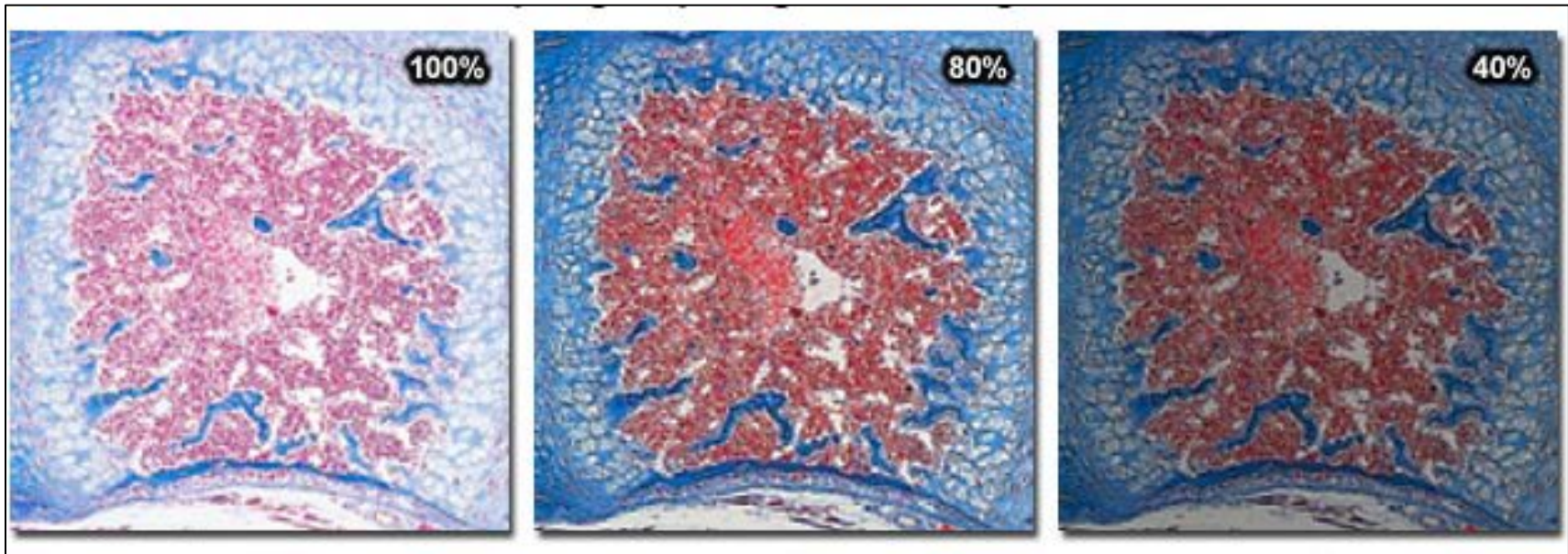



We make it visible.

<http://zeiss-campus.magnet.fsu.edu/index.html>